


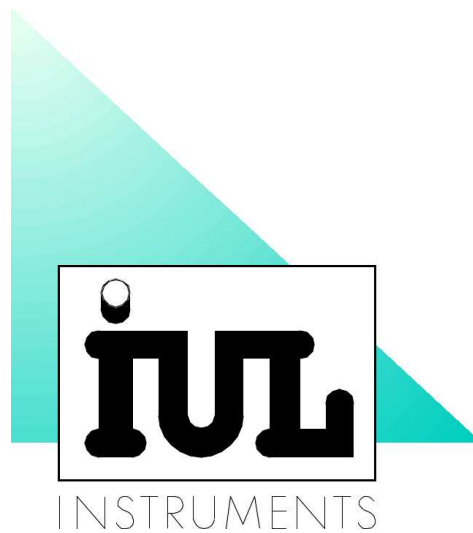

MASTICATOR COMPACT CLASSIC


USER'S GUIDE

DOC. No. 156R01

HOMOGENEIZER


MASTICATOR
COMPACT
CLASSIC

INSTRUCTION
MANUAL

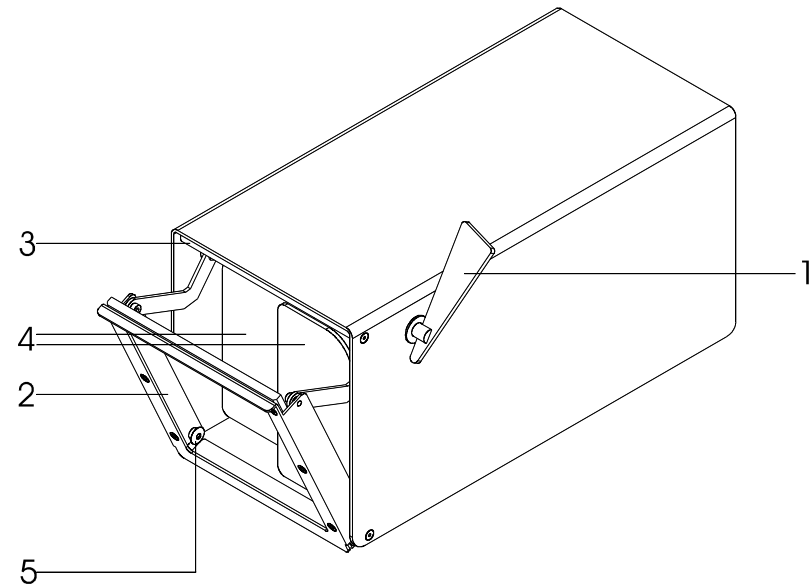


FRONT PART

- 1. Opening and closing lever.
- 2. Sample compartment lid.
- 3. Bag closing bar.
- 4. Homogenization blades.
- 5. Lid bracket.

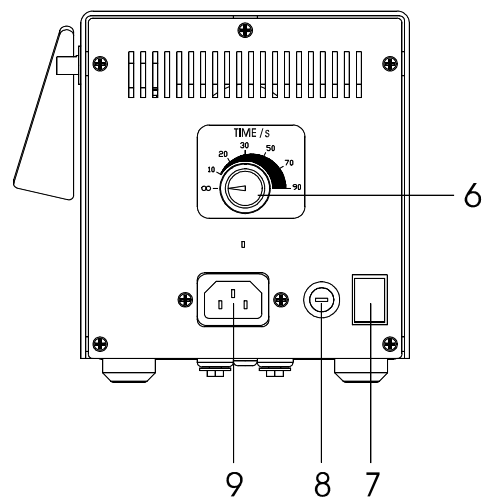

REAR PART

- 6. Time switch
- 7. Main power switch ON/OFF.
- 8. Fuse.
- 9. Mains lead connection.

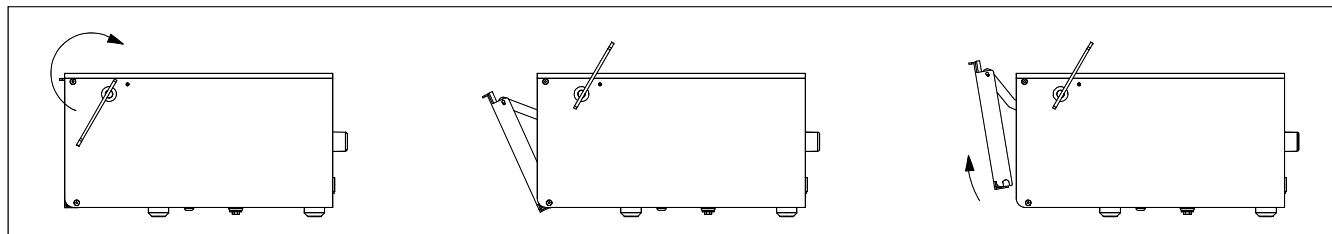

INDEX

1. INTRODUCTION	1
2. TECHNICAL SPECIFICATIONS	2
3. INSTALLING AND STARTING THE INSTRUMENT.....	3
4. OPERATION AND USE THE CONTROLS	4
5. BEFORE NOTIFYING THE TECHNICALSERVICE	5
6. ACCESSORIES AND OTHER RELATED INSTRUMENTS	6
7. WARRANTY CONDITIONS.....	7

USE OF THE MANUAL

All the letters, which are circled; (A), (B), (C) etc. refer to the stages, which must be carried out by the operator. All the numbers; (no.1), (no.2), (no.3), etc. refer to the main illustration on the back cover.

IMPORTANT

Read this little manual carefully before starting the apparatus. It will not take you very long and it will allow you to get to it more quickly and more conveniently.

.

1. INTRODUCTION

In many microbiological analysis processes it is essential to have available a liquid sample which has a bacterial content which is representative of a solid, semi-solid, paste and pulverulent, etc. original material. The liquid sample is easy to cultivate or to process microbiologically whereas the original material would not allow reliable results to be obtained.

In general the homogenization process involves the original substance and the support liquid or diluent being brought into close contact, in order to establish a balance in the microbiological content of both of these.

This process facilitates the use of reusable material, which must be washed and sterilized so as not to invalidate consecutive samples.

WHAT IS THE MASTICATOR?

The Masticator is a simple and yet robust instrument, designed essentially to work for prolonged periods without maintenance, in any kind of the environment. Initially devised for the homogenization of different materials in order to take representative samples from them, it is currently used for other applications related to microbiological analysis.

In a short period of time, and with the minimum amount of noise, and without the necessity of cleaning between samples, the Masticator carries out a profound homogenization which makes it possible to obtain a diluent liquid which is perfectly representative of the sample which is to be analyzed.

Its working principle lies in the use of bags made of strong polyethylene or other plastics, inside which the original sample and the diluent liquid are placed. The Masticator acts upon the sample and the diluent through the sides of the bag, which, at this moment, is perfectly sealed, without either soiling or contaminating any part of it whatsoever. The action of the blades compresses and allows the periodic expansion of the sample and the diluent for a pre-determined period of time. This operating time is controlled by a time switch, which can be adjusted by the operator. In spite of this, a process which is in progress can be interrupted by activating the lever in order to open the homogenization compartment.

The parts, which are in contact with the sides of the bag, are made of stainless steel, and both the stopper and the homogenization blades are perfectly accessible in order to allow them to be eventually cleaned.

2. TECHNICAL SPECIFICATIONS

- Maximum capacity per bag: 100 ml. (sample + diluent)
- Strokes/s: 10
- Time switch minimum: 10 seconds
- Time switch maximum: 90 seconds
- Position of uninterrupted operation: Yes
- Power supply: 230V, 0.3 A, 50/60 Hz ,CAT. No. 85
115V, 0.6 A, 60 Hz ,CAT. No. 87
100V, 0.6 A, 50/60, Hz CAT. No. 86
- Dimensions: 150 x 150 x 300 mm.
- Weight: 7 Kg.

3. INSTALLING AND STARTING THE INSTRUMENT

The Masticator requires no special conditions for its installation. Its robustness and its reduced size allow it to be located in any part of the laboratory whatsoever. The only thing which it is necessary to take into account is the proximity of a power source and a surface large enough for it to be correctly supported.

N.B. DO NOT EXPOSE THE SAMPLE COMPARTMENT TO THE AIR. IF THE TABLE IS LARGE ENOUGH, IT WILL NEVER BE POSSIBLE, NOT EVEN BY ACCIDENT, FOR A HAND TO BE PASSED BENEATH THE LID OF THE HOMOGENIZATION COMPARTMENT WHEN THE MASTICATOR IS IN OPERATION, IN THIS WAY AVOIDING THE POSSIBLE TRAPPING OF THE FINGERS.

The Masticator should be received in perfect condition, it was been tested at the factory and the packaging means that it can stand up to the roughest treatment.

Once the instrument has been placed upon the work table, plug the mains lead into its socket, 9, and, at the base, locate the time switch in its minimum position (ten seconds) and switch on the main power switch ON/OFF, 7. The Masticator will start up, stopping after approximately ten seconds.

By means of the lever, 1, open the homogenization compartment.

Remove the lid of the compartment, do this by means of a light movement pulling it upwards and detaching it from its interior brackets 5.

Now look at the inside, clean off any possible remains of the packaging and make sure that the homogenization compartment is open in its lower part.

This design permits the easy elimination of any remains that could be left over after the accidental breaking of a bag on the inside of the homogenization compartment.

Tilt the lid to its working position and fit the lower angle onto its brackets. The instrument is now ready for operation

4. OPERATION AND USE THE CONTROLS

- Ⓐ Set the time switch to halfway along its possible setting. Thirty seconds is an enough time for most samples.
- Ⓑ Introduce the sample and the solvent into a homogenization bag. Do not exceed 100 ml in total.
- Ⓒ Place the bag into the compartment, the lower edge where the lid, 2, is held down by the brackets 5, is an end where the bag may be rested.
- Ⓓ With the mouth of the bag very flat, held with one hand and supported by the rim of the lid, let out 3 cm. and go on to close the compartment with the lever, 1. A moment before the engine starts up the bag will have been hermetically sealed by the pressure of the lid upon the upper bar, 3. The mouth of the bag is trapped between the bar and the lid and is only free when the lever, 1, is activated again.
- Ⓔ Once the period of the pre-determined time with the time switch is up, the motor will stop. Hold the flap of the bag and open the compartment with the lever, 1. Take out the bag.

5. BEFORE NOTIFYING THE TECHNICAL SERVICE

If the Masticator does not start up check the connection of the mains lead and the fuse, (no.8), situated on its rear side. If both of these are correct contact the technical service in your country or:

IUL S.A.

Torrent de l'Estadella 22

08030 BARCELONA - SPAIN

Phone: +34 -932740232


Fax: +34 - 932740144

PLEASE NOTE:

THIS INSTRUMENT WILL NEED TO BE KEPT CLEAN BY THE USER. USE A CLOTH MOISTENED WITH WATER BUT NEVER FORGET TO DISCONNECT THE INSTRUMENT FROM THE POWER SUPPLY BEFORE BEGINNING TO CLEAN IT. THE STAINLESS STEEL COMPARTMENT CAN BE CLEANED WITH DISSOLVENTS OR OTHER MORE ACTIVE AGENTS, WITHOUT ANY DANGER OF ITS BEING DAMAGED IN ANY WAY

6. ACCESSORIES AND OTHER RELATED INSTRUMENTS

The **MASTICATOR** is designed not to need any kind of periodic maintenance.

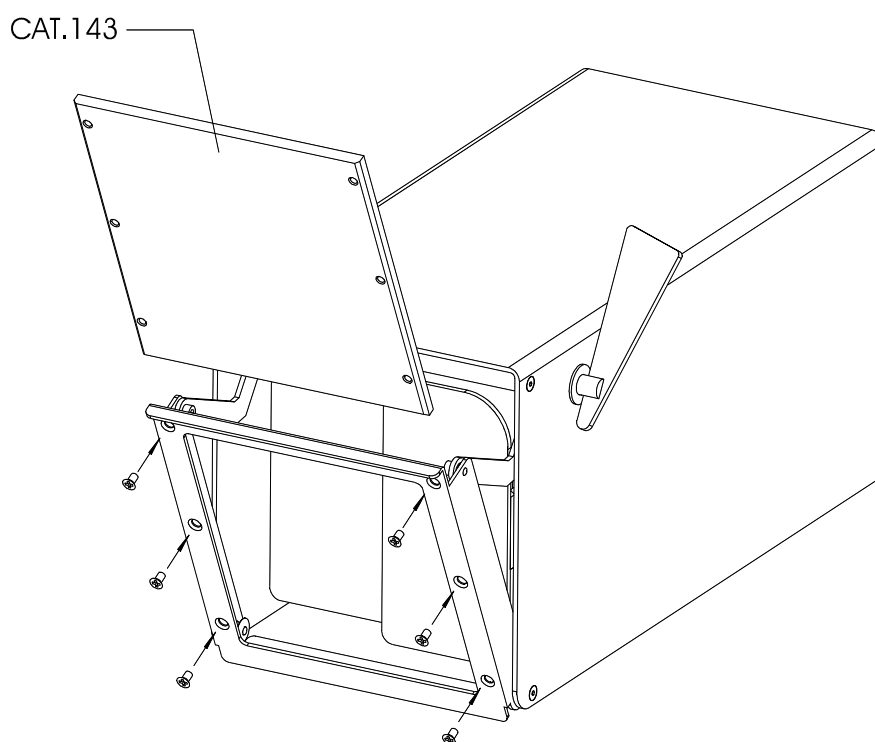

Replacing the polycarbonate lid

CAT. 143 Polycarbonate lid

The **DELTA DILUTOR** and the **PORTA-BAG** are also included in this section.

The **DELTA DILUTOR** is a diluter especially designed for carrying out the dilution of microbiological samples directly into sterile homogenizer bags. The system automatically prepares your desired dilution based upon sample weight, liquid flow rate and predefined dilution factor.

The **PORTA-BAG** is a device which allows six bags filled with the sample and the diluent to be held in a vertical position without the fear of spilling the contents. The tray which holds the pincers and these themselves are made of stainless steel.

7. WARRANTY CONDITIONS

LIMITED WARRANTY

IUL S.A. warrants the Instrument, disposables and replacement items, to the original "Purchaser" for a period of two (2) years after date of receipt by the end user, against defects in material and workmanship. **IUL S.A.** agrees to correct either by repair or, at its election, by replacement, any such defect found on examination to have occurred, under normal use and service, during such one (1) year period provided IUL is promptly notified in writing upon discovery of such defect to the address of one of its workshops noted in the User's Guide.

Upon discovery of such defect, the instrument should be sent to an IUL workshop noted in this User's Guide, postage prepaid.

EXCLUSIONS

IUL S.A. shall not be liable as to any defect arising from abuse of the Instrument, failure to operate and maintain the Instrument in accordance with the User's Guide supplied with the Instrument, repair, service, alteration or modification of the Instrument by any person other than the authorized service personnel of IUL, modification, change or reuse of the disposables supplied by IUL for use in the Instruments, operation outside of the environmental specifications for the product listed in the User's Guide or accidental damage to the Instrument.

NO OTHER WARRANTIES

The warranty set forth above is exclusive and no other warranty, whether written or oral, is expressed or implied. IUL and its suppliers specifically disclaim any implied warranties including merchantability and fitness for a particular purpose.

NO LIABILITY FOR CONSEQUENTIAL DAMAGES

IUL or its suppliers shall not be liable as to any damage whatsoever (including without limitation, special, incidental, consequential, or indirect damages for personal injury, loss of business profits, business interruption, loss of business information, or any other pecuniary loss) or defect arising from abuse of the Instruments, failure to operate and maintain the Instruments in accordance with the User's Guide, repair, service, alteration or modification of the Instruments by any person other than authorized service personnel of IUL, or modification, change or reuse of the disposables supplied by IUL for use in the Instruments

IUL S.A
Torrent de l'Estadella 22
08030 Barcelona-Spain
Tel. (34) 932740232 – Fax. (34) 932740144

www.iul-inst.com